


Figure S11

Blm

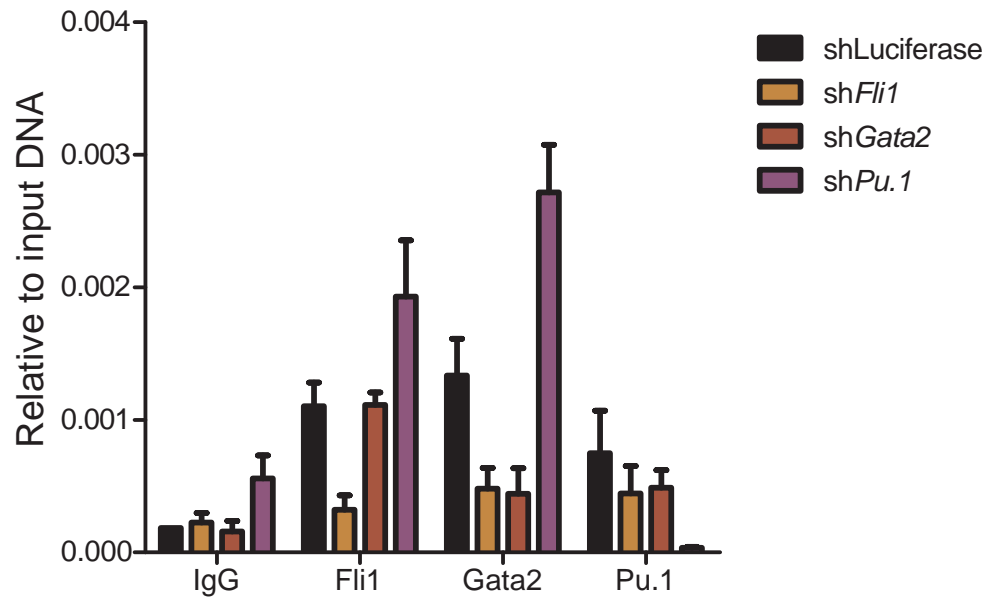

Tspan2

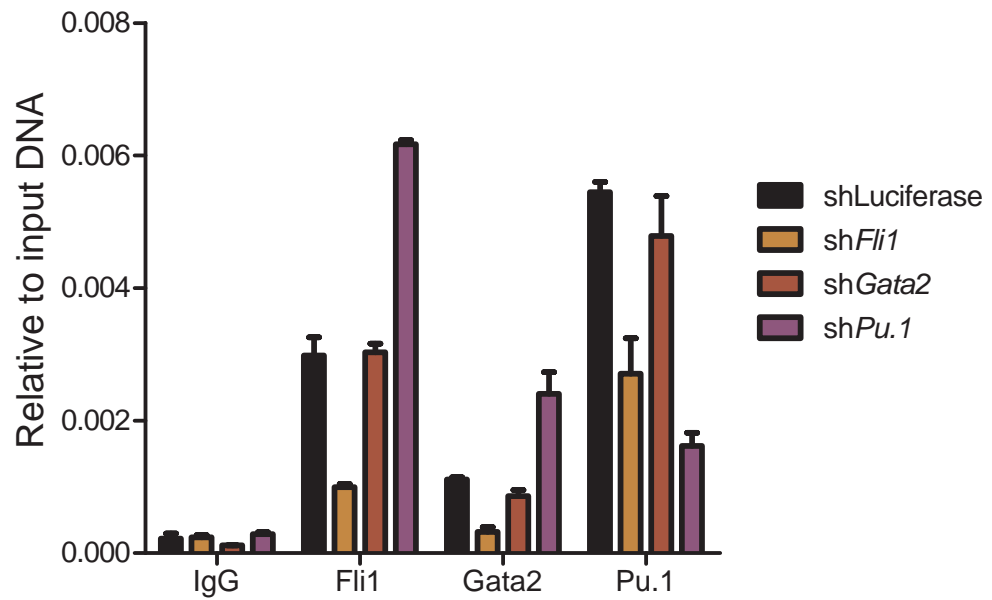


Figure S11 – Real-time PCR analysis of ChIP for FLI1, GATA2 and PU.1 following knock-down of the same genes in regions associated to the genes *Blm* and *Tspan2*. Experiments were performed in MST mast cell line and results are expressed relative to total DNA input. Primers used for amplifications were: *Blm* (5' GGGATGAGTCACTGCCAACT 3' and 5' ATTTGCCTTGTTTCCTGCTG 3'); *Tspan2* (5' GCCCCATCTAGGTCTGTTCA 3' and 5' ACTCTGGCATTTCAGGCACTT 3').